



TOSOH

IVF SOLUTIONS



**E2 Results in 18
to 20 Minutes**

Why Wait?

*Tosoh Immunoassay Analyzers
Deliver Accurate Results Quickly*

Fast - First Result in 18 to 20 Minutes

Efficient - Simple Automated Testing


Accurate - Correct Results the First Time

Consistent - Common Unit Dose Test Cup Reagent

Cost Saving - 90 Day Calibration Stability for all Fertility Assays

Comprehensive - Full Immunoassay Menu Available, Including a Full Menu of Reproductive Hormone Assays

Biotin Free - All Tosoh's AIA-PACK® Immunoassays are Free from Biotin Interference

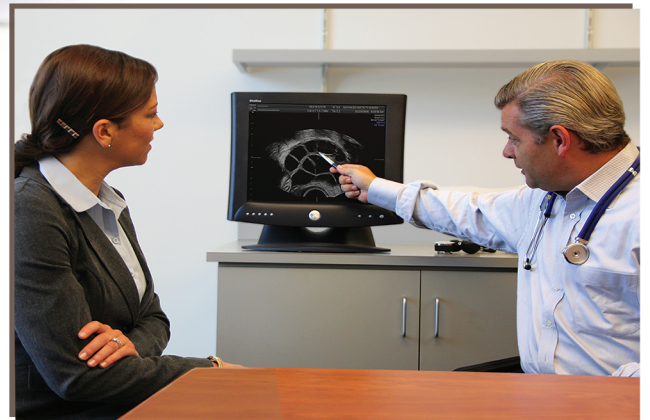

*Tosoh's fast results allow review of ultra-sound
and E2 on the same visit.*

TOSOH BIOSCIENCE

Unit Dose Test Cup Technology

- Single, unitized cups use a dry reagent format that requires no pre-mixing, no pre-measuring, no on-board refrigeration and no waste
- Dry reagent format ensures calibration stability of up to 90 days for most assays
- Reagents have a one year shelf life from date of manufacture
- The AIA test cups and trays are labeled with code and lot number for automated scheduling and inventory management on AIA-2000 and AIA-900 sorter models
- ST AIA-PACK reagents provide a faster turnaround time due to a 10 minute incubation protocol
- Extensive test menu provides a wide range of immunoassay tests including thyroid, cardiac markers, tumor markers, reproductive hormones, anemia markers, metabolic markers, diabetes, kidney markers and more
- Biotin does not affect Tosoh assays because the methodology does not utilize the biotin-streptavidin pair to isolate the signal



Available for Labs and Workloads of All Sizes



AIA-360
36 Tests Per Hour



AIA-900 Benchtop
90 Tests Per Hour



AIA-900 Loader Model
90 Tests Per Hour



AIA-900 with 9 tray sorter
90 Tests Per Hour
Also available with 19 Tray Sorter



AIA-2000
200 Tests Per Hour

Choose Wisely. Choose Tosoh.



TOSOH BIOSCIENCE

Tosoh Bioscience, Inc.
3600 Gantz Rd
Grove City, OH 43123
Tel: (800) 248-6764
Fax: (650) 615-0415
www.tosohbioscience.us

"AIA-PACK" and "AIA" are registered trademarks of Tosoh Corporation